
Mepilex[®] Border Lite

WITH SAFETAC[®] TECHNOLOGY

STERILE EO

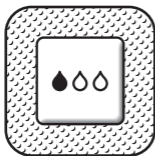Single use
Usage unique

CE 0086

MÖLNLYCKE
HEALTH CARE

USA 1-800-882-4582
Canada 1-800-494-5134

Manufacturer

Mölnlycke Health Care AB

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Mepilex® Border Lite with Safetac® Technology

Self-adherent soft silicone foam dressing

Product description

Mepilex Border Lite consists of:

1. a soft silicone wound contact layer [Safetac®] on a polyurethane film carrier
2. an absorbent pad in two layers: a polyurethane foam and a viscose/polyester non-woven spreading layer

3. an outer polyurethane film which is vapour permeable and waterproof

Safetac is a unique and a patented adhesive technology that minimises pain to patients and trauma to wounds and the surrounding skin.

Mepilex Border Lite is a thin, highly conformable self-adherent dressing that absorbs exudate, maintains a moist wound environment and minimises the risk for maceration.

As Mepilex Border Lite maintains a moist wound environment, supporting debridement, there might be an initial increase in the wound size. This is normal and to be expected.

Intended use

Mepilex Border Lite is designed for the management of a wide range of non/low exuding wounds, such as leg and foot ulcers, pressure ulcers and traumatic wounds, e.g. abrasions, finger injuries, blisters and skin tears.

Mepilex Border Lite can also be used as protection of compromised and/or fragile skin.

Instructions for use

1. Cleanse the wound in accordance with normal procedures.
2. Dry the surrounding skin thoroughly.
3. Remove the release films and apply the adherent side to the wound.
Do not stretch.
4. For best results, Mepilex Border Lite should overlap the dry surrounding skin by at least 1-2 cm for the smaller sizes (sizes up to 12.5x12.5 cm) and 5 cm for the larger sizes in order to protect the surrounding skin from maceration and fixate the dressing securely.

Mepilex Border Lite may be left in place for several days depending on the condition of the wound and surrounding skin, or as indicated by accepted clinical practice.

A change in dressing regimen can result in an initial increased level of exudates, which temporarily may require an increased change frequency.

Mepilex Border Lite can be used under compression bandaging

Mepilex Border Lite can be used in combination with gels

Precautions

Do not use on patients with known sensitivity to the dressing or its components.

In the case of signs of clinical infection, consult a health care professional for adequate infection treatment.

Do not use Mepilex Border Lite together with oxidising agents such as hypochlorite solutions or hydrogen peroxide.

Storage and disposal

Mepilex Border Lite should be stored in dry conditions below 35°C (95°F).

Disposal should be handled according to local environmental procedures.

Other information

The polyurethane foam used in the product may change colour to more yellow when it is exposed to light, air and/or heat. The colour change has no influence on product properties when used before expiry date.

Do not reuse. If reused performance of the product may deteriorate, cross contamination may occur.

Sterile. Do not use if inner package is damaged or opened prior to use. Do not re-sterilise.

If the product is used after expiry date product properties cannot be ensured.

Mepilex and Safetac are registered trademarks of Mölnlycke Health Care AB.

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CANADA Mölnlycke Health Care Inc.,
2010 Winston Park Drive, Suite 100,
Oakville, Ontario, L6H 5R7
Tel: +1 905 829 1502

CHINA Mölnlycke Health Care
(Shanghai) Co., Ltd.
RM 629, No.8 Hua Jing Road,
Wai Gao Qiao FTZ Shanghai, China
Tel: +86 10 5128 8571

CZECH REPUBLIC Mölnlycke Health Care s.r.o.
Hajkova 2747/22, 130 00 Prague 3
Tel: +420 221 890 517 **[Reception]**
Tel: +420 221 890 511 **[Customer Service]**

DENMARK Mölnlycke Health Care ApS,
Gydevang 33, DK-3450 Allerød
Tel: +45 48 16 82 68
Customer Service Tel: +80 886 810

ESTONIA Mölnlycke Health Care OÜ,
Vabaõhumuuseumitee 4-67, EE-13522 Tallin
Tel: +372 671 1520

FINLAND Mölnlycke Health Care Oy,
Pitäjänmäentie 14, 3 krs, FI-00380 Helsinki
Tel: +358 201 622 300
Customer Service Tel: 0800 113 442

FRANCE Mölnlycke Health Care SAS,
13 Allée du Château Blanc
59290 WASQUEHAL
Tel: +33 320 122 555
Customer Service Tel: 0800 910 292

GERMANY Mölnlycke Health Care GmbH
Grafenberger Allee 297
DE-40237 Düsseldorf
Tel: +49 211 92 08 80
Customer Service Tel: 0800 186 21 80

HUNGARY Mölnlycke Health Care Kft
HU-1134 Budapest, Dévai u. 26-28. III. em.
Tel: +36 1 477 3080

INDIA Mölnlycke Health Care India Pvt Ltd
Universal Business Centre
Second Floor, Universal Trade Tower
Sector 49, Gurgaon Sohna Road
Gurgaon 122 018, Haryana
Tel: +91 124 4696 008

ITALY Mölnlycke Health Care s.r.l.,
Via Marsala 40/C, IT-21013 Gallarate (VA)
Tel: +39 0331 714 011

JAPAN Mölnlycke Health Care KK,
4F Concieria Nishi-Shinjuku Tower's West
6-20-7 Nishi-Shinjuku,
Shinjuku-Ku, Tokyo 160-0023
Tel: +81 3 6914 5004

KOREA Mölnlycke Health Care Korea Co., Ltd
2nd floor, Hansol Richville.
46 Wiryeseong-daero,
Songpa-gu, Seoul, 05627
South Korea
Tel: +82 2 3789 1402

LATVIA Mölnlycke Health Care AB,
Dzelzavas 120 M, Riga LV-1021
Tel: +372 6711 522

LITHUANIA Mölnlycke Health Care AB,
J. Galvydzio g. 3, LT 08236 Vilnius
Tel: +370 5 2274 58 37

NETHERLANDS Mölnlycke Health Care B.V.
Postbus 3196, 4800 DD Breda
Tel: +31 76 521 9663

NORWAY Mölnlycke Health Care AS
Postboks 6229 Etterstad, NO-0603 Oslo
Tel: +47 22 70 63 70
Customer Service Tel: 800 161 78

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Mölnlycke Health Care AG
Mölnlycke Health Care Ltd.
USA Mölnlycke Health Care US, LLC